



When to select <99% Filtering Efficiency on a Condensation Nuclei Counter (CNC) instrument during Respirator set up

A CNC instrument grows the microscopic particles then, counts the particles inside a respirator and compares them to the particles counted on the outside of the respirator (C_{out}/C_{in}), producing a fit factor. The CNC instrument counts all particles (in the specified size range) and operates on the premise that all particles within the respirator sample have come via the face seal leakage pathway.

When setting up a respirator in a CNC instrument it is critically important to know when to select the <99% filtering efficiency option or not. Whether it needs to be selected is based upon the filtering efficiency of the particulate filter media used in the fit test. Refer to Appendix A for screenshots of different CNC instrument software and screens.

If the filtering efficiency of the particulate filter media being used in the fit test is $\geq 99\%$ such as P3, N100, P100, R100, FFP3 or P3 R (but not limited to), the <99% filtering efficiency DOES NOT need to be selected.

If the filtering efficiency of the particulate filter being used in the fit test is $< 99\%$ such as P1, P2, N95, R95, P95, KN95 FFP1 or FFP2 (but not limited to), the <99% filtering efficiency DOES need to be selected.

These particulate filter standards (AS/NZS, NIOSH, EN etc) ratings allow a certain amount of the standard test aerosol (being the most penetrating particle size) to penetrate the filter and still meet certain ratings. When <99% filtering efficiency is selected the CNC instrument counts a specific size of particles (reduced range) rather than a broader range of sizes (which is why the particle count goes down when this is selected). When the <99% filtering efficiency button is selected, a reduced range of particles is counted by the CNC instrument, eliminating the most penetrating particle sizes (MPPS). The particle sizes that are counted are known to be captured at >99% for these rated particulate filter levels. As a result, we are truly testing the seal of the respirator as a gap in the seal is the only way these particles could have entered and been counted in the facepiece sample.

If the <99% filtering efficiency is not selected when fit testing using P1, P2, N95, R95, P95, KN95, FFP1 or FFP2 filter during the fit test, the MPPS particles will be penetrating/passing through the filter (which is allowable per the standard rating). This fit test method assumes that all counted particles have gone through the face seal pathway only, not through the filter. This will lead to a false negative result (i.e., potentially failing good fitting respirators and not knowing it).



A gas/vapour cartridge cannot be used for a fit test UNLESS a suitable rated particulate filter has been attached/stacked to it. For a fit test that is only conducted using a gas/vapour cartridge, the result should be voided/discarded and another fit test with an appropriate particulate filter should be conducted.

Please refer to the table below to assist in when to select the < 99% filtering efficiency option in your CNC instrument respirator set-up.

Respirator Style	Common Particulate Filter Ratings Used in Australia (but not limited to)	Select <99% Filtering Efficiency CNC Machine
Disposable/Filtering Facepiece	P1, P2 - AS/NZS 1716 N95, P95, R95 – NIOSH US FFP1, FFP2 – European Standards KN95 – Chinese Standard	Yes
	Note: Under AS/NZS 1716 it does not have a >99% disposable/filtering facepiece rating N99, N100, R99, R100, P99, P100 – NIOSH US FFP3 – European Standards	No
Reusable/Elastomeric Half Face*	P1, P2 - AS/NZS 1716 N95, P95, R95 – NIOSH US P1 R, P2 R – European Standards	Yes
	P3 - AS/NZS 1716 N99, N100, R99, R100, P99, P100 – NIOSH US P3 R – European Standards	No
Reusable/Elastomeric Full Face	Note: To fit test a full face, a filter of $\geq 99\%$ must be used. A <99% filter cannot be used e.g. P2, N95, KN95 etc. P3 - AS/NZS 1716 N99, N100, P99, P100, R99, R100 – NIOSH US P3 R – European Standards	No

*It is common practice that many fit testers will use a $\geq 99\%$ particulate filter e.g. P3, on Reusable/Elastomeric Half Face respirators as generally, they will also connect to the same brand/connection type reusable/elastomeric full face as part of their fit testing equipment. If fit testing a Reusable/Elastomeric Half Face and changing between <99% and $\geq 99\%$ particulate filters, be sure to review your respirator set-up each time is correct to the filter efficiency being used.

Should you have any questions, please email respfit@aioh.org.au


Version: v2 Document Number: AIOH-RESPFIT-17

Next Review Date: September 2026




Appendix A – Pictures of Different CNC Software and CNC Machine Screens where to select <99% filtering efficiency

TSI FitPro Ultra

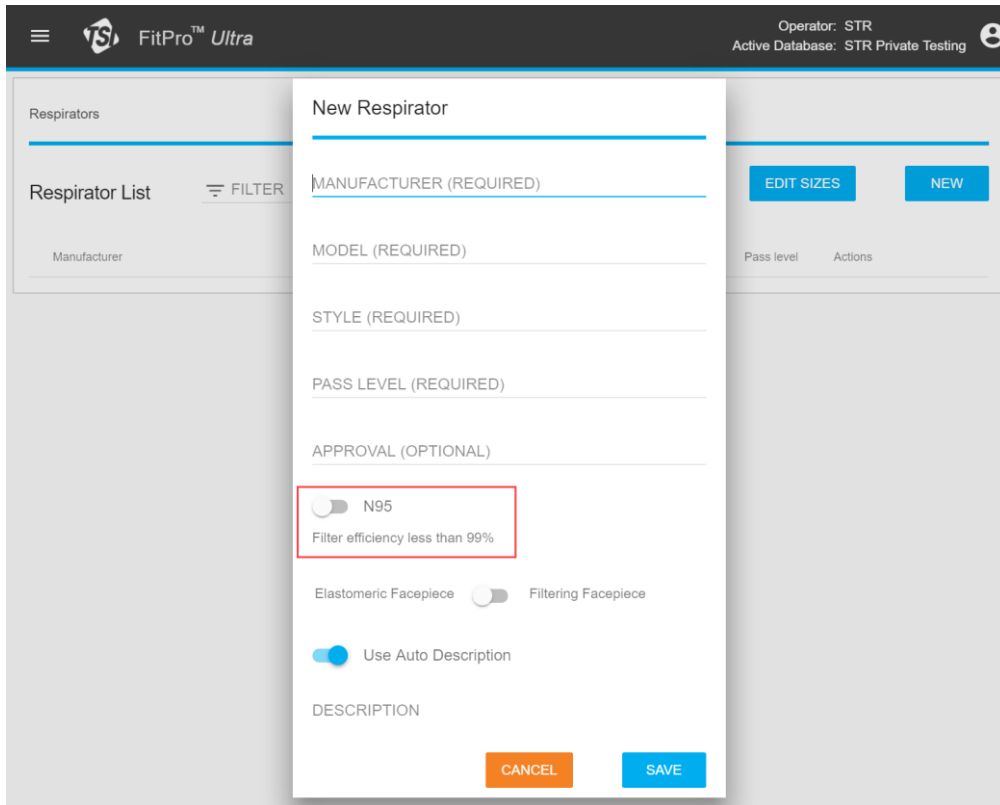

TSI FitPro+

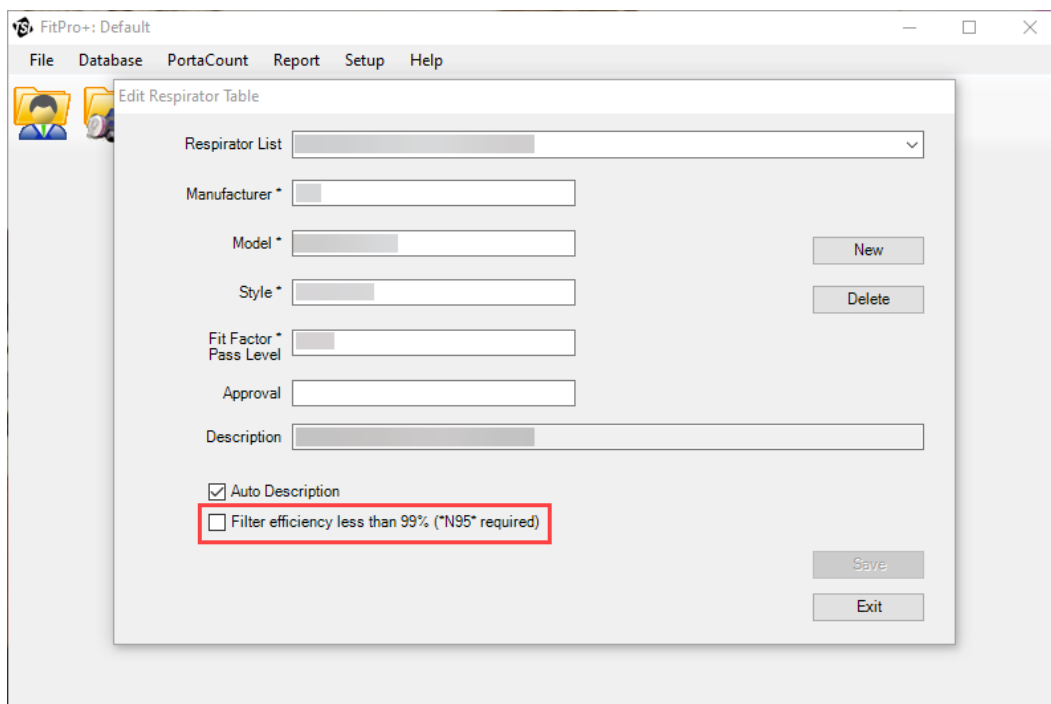



TSI PortaCount Respirator Fit Tester 8038 Machine Screen

The screenshot shows a software interface titled "Create New Entry" for a "Respirator List". The form includes several input fields: "Manufacturer *", "Model *", "Style *", "Pass Value *" (set to 100), and "Approval". The "Approval" field has a dropdown menu with "N95" selected and highlighted by a red box. A "Save" button is located to the right of the "Manufacturer" field, and a "Cancel" button is at the bottom right. A note at the bottom left states "* Required Fields".

AccuTec-IHS AccuFIT DataManager

The screenshot shows a web-based interface titled "RESPIRATOR DATABASE PORTAL - NEW ENTRY". It features a large "SAVE" button at the top. Below it is a table with four columns: "Manufacturer", "Model", "Type", and "Pass Value". At the bottom, there is a checkbox labeled "Filter efficiency less than 99%" which is highlighted by a red box. An "EXIT" button is located at the bottom right.



AccuTec-IHS AccuFIT9000 Pro Machine Screen

Respirators | Demo 09/04/2020 10:24

RESPIRATOR DATABASE PORTAL - NEW MASK

Manufacturer	Model	Type	Pass Value
<input type="text"/>	<input type="text"/>	<input type="text"/>	<input type="text"/>

Filter efficiency less than 99%